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**SPORT RELAX LB8787**

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**28130 Saint-Piat**

**France**

**INSTRUCTION MANUAL LB8787 DEEP TISSUE MASSAGER**  
**SPORT RELAX**

## IMPORTANT SAFETY INSTRUCTIONS

Please read these instructions thoroughly before first use, even if you are familiar with this type of product. The safety instructions are intended to reduce the risk of fire, electric shock and injury when correctly adhered to. We recommend keeping this instruction manual in a safe place for future reference. If you pass the appliance on to anyone else, please include this instruction manual.

### TECHNICAL SPECIFICATIONS

Model: LB8787

NON-Waterproof

5V  $\equiv$  2A

Battery Li-ion 1200 mAh DC 7.4V

Power: 30W

Motor: 1800 – 3200 RPM

Type-C charging cable

Charge time: 200 minutes

Cordless use time: 90 minutes

Class III product

### DESCRIPTION

This appliance is exclusively intended for massage treatment on the legs, shoulders, arms and back region (excluding the spine). The Deep Tissue Massager is NOT intended for use on the head or face. DO NOT use it on pets.

The Deep Tissue Massager is NOT a medical device but a massage treatment ONLY.

This appliance is ONLY suitable for personal use, NOT for medical or commercial purposes.

The Deep Tissue Massager may ONLY be used for the purposes for which it is designed and, in the manner, specified in the instructions for use. Any other use might lead to damage of the appliance, property or personal injury. The manufacturer will accept no liability for damage or injuries resulting from improper use or non-compliance with this manual.

## GENERAL SAFETY PRECAUTIONS

In order to reduce the risk of fire, electrical shock or injury to persons or property, follow these basic safety precautions when using any electrical product:

Before using the appliance for the first time, ensure all packaging and stickers are removed. Dispose of them properly, keeping out of reach of children. DO NOT remove the rating label from this product. It contains important safety information.

**WARNING! RISK OF ELECTRIC SHOCK!** DO NOT immerse the appliance, USB charging cable or adapter into water or any other liquid. If the appliance falls into water, unplug it before touching it. This appliance is NOT suitable for a use in the bath or shower.

**WARNING!** When the appliance is used in a bathroom, switch it off after use since the proximity to water presents a hazard even when the appliance is switched off.

**WARNING!** This device is NOT suitable for cleaning under an open water tap.

**WARNING!** The device must ALWAYS be disconnected from the USB charging cable before cleaning.

When not in use, ALWAYS switch off the appliance. Do the same before cleaning the appliance or attaching or removing any attachments.

This appliance can be used by children aged from 8 years and above, and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge ONLY if they have been given instruction concerning the correct use of this appliance or are using it while supervised and understand the hazards involved.

This appliance is NOT a toy. Children should NOT play with the appliance.

Cleaning of this appliance should NOT be carried out by children.

**WARNING!** NEVER leave the appliance unattended while in use or plugged in.

Prior to charging the device, check that your power source has the same voltage, frequency and rating as indicated on the appliance.

Ensure that your hands are dry when handling the appliance. DO NOT operate the appliance with wet hands.

DO NOT let the USB charging cable to hang over the edge of a table or counter, touch hot surfaces or become knotted. Protect it from sharp edges.

ALWAYS place the appliance on a solid, flat and dry surface when operating.

**WARNING!** DO NOT operate the appliance when it is covered by a blanket or similar. This can lead to overheating and risk of fire.

DO NOT use the appliance with a timer or remote control.

Inspect the appliance and USB charging cable regularly. DO NOT use the device after it malfunctions, has been dropped or appears damaged in any way. Return the complete product to the place of purchase for inspection, replacement or refund.

**WARNING!** Should this product malfunction, DO NOT attempt to repair it. This product has no user-serviceable parts or components.

DO NOT open or tamper with the appliance casing (housing) under any circumstances.

Ensure that all attachments and accessories are properly fixed, to prevent them from falling into the water.

DO NOT use the appliance with any accessories not recommended or supplied by the manufacturer. This will invalidate your warranty and may damage the appliance or injure the user.

DO NOT use this appliance for any other purpose than indicated. This appliance is intended solely for use in private households.

Commercial use invalidates the warranty. The supplier cannot be held responsible for injury or damage caused when using the appliance for any other purpose than intended.

DO NOT use this appliance outdoors. It is intended for indoor use ONLY.

## BATTERY

This appliance contains batteries that are non-replaceable.

**WARNING!** DO NOT dismantle, open or shred the battery.

**WARNING! DANGER OF EXPLOSION!** Protect the battery against heat, e.g, against continuous intense sunlight or fire, and against water and moisture! Also protect the battery from cold!

## USB CHARGING CABLE

**WARNING!** For the purposes of recharging the battery, ONLY use the detachable USB charging cable provided with this appliance.

DO NOT attempt to charge the appliance with any other USB charging cable than the one supplied.

The appliance and USB charging cable supplied with it are specifically designed to work together.

ONLY use the provided USB charging cable or the one recommended by the manufacturer.

If the USB charging cable is damaged, it shall be scrapped. Disconnect the USB charging cable from the power supply before cleaning, maintenance, storage and transportation.

DO NOT touch the USB charging cable and adapter with wet hands.

DO NOT let the USB charging cable hang over sharp edges, DO NOT squeeze or bend it. Keep it away from hot surfaces.

**WARNING! RISK OF ELECTRIC SHOCK!** Protect the USB charging cable against moisture and wet conditions. Charge the battery in an ambient temperature of 18 - 24 °C. DO NOT charge the battery in temperatures below 0 °C and above 40 °C. Charge the battery indoors ONLY.

DO NOT place any object on top of the USB charging cable, as it could cause overheating.

### IMPORTANT MEDICAL WARNINGS

Under NO circumstances use the device for massaging the area around the heart.

NEVER use the device on swollen, burnt, inflamed or injured skin and body parts.

DO NOT use the device if you suffer from a medical abnormality or injury in the area to be massaged (e.g slipped disc, open wounds).

DO NOT use the device on cuts, capillaries, varicose veins, acne, couperose, herpes or other skin diseases.

DO NOT use the device during pregnancy.

DO NOT use the device on the head.

DO NOT use the device on the face (eyes), the larynx or other particularly sensitive body parts.

DO NOT use the device on the spine.

DO NOT use the device after the consumption of medicines or alcohol (reduced cognition!).

**You should consult your doctor before using the device, especially:**

If you suffer from a severe illness.

If you have a pacemaker, implant or other aids.

If you suffer from thrombosis.

If you suffer from diabetes.

If you are suffering from pain, the cause of which is unknown.

## SAFETY PRECAUTIONS FOR YOUR LITTLE BALANCE DEEP TISSUE MASSAGER

Use the Deep Tissue Massager ONLY for purposes described in this instruction manual.

**WARNING! - RISK OF ELECTRIC SHOCK!** DO NOT immerse the appliance, USB charging cable or adapter into water or any other liquid. If the appliance falls into water, unplug it before touching it.

**WARNING!** This appliance is NOT waterproof. It is NOT suitable for a use in the bath or shower.

**WARNING!** When the massager is used in a bathroom, switch it off after use since the proximity to water presents a hazard even when the massager is switched off.

DO NOT place or store the massager where it can fall into water or any other liquid.

**WARNING!** This device is NOT suitable for cleaning under an open water tap.

Please ensure that your hands are dry when handling your massager. DO NOT operate the appliance with wet hands.

**WARNING!** NEVER leave the appliance unattended while in use or plugged in.

ALWAYS switch off the massager after use or before cleaning it.

Regularly inspect the appliance and USB charging cable for any damage. If found damaged in any way, immediately cease use of the appliance.

NEVER use the massager after it malfunctions or if it has been dropped or damaged in any way.


If you do find any damaged items, then DO NOT use the massager. Contact your place of purchase for advice.

NEVER disassemble the massager.

## PRODUCT OVERVIEW

## 4 REPLACEABLE MASSAGE HEADS

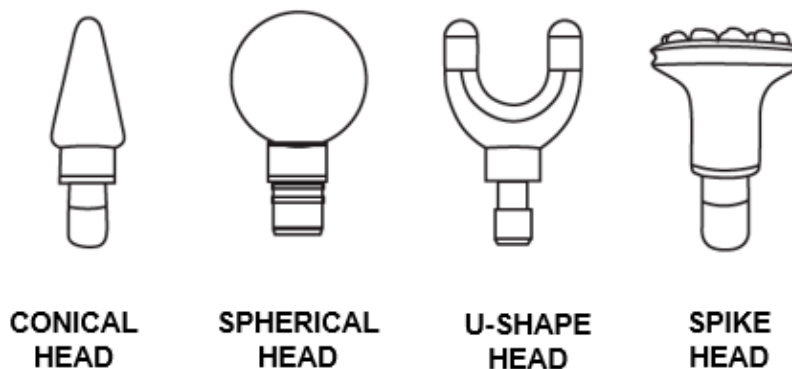**CONICAL HEAD:** Best suited for treating deep muscle tissue and for targeted spot massages.

**SPHERICAL HEAD:** Best suited for treating large areas and large muscle groups.

**U-SHAPED HEAD:** Best suited for treating deep muscle tissue in the lower back area (excluding the spine) or the Archilles tendon.

**SPIKE HEAD :** Best suited for treating large areas on the entire body (legs and deep stimulation of the buttocks).

## BEFORE FIRST USE


**WARNING!** Packaging materials are not a toy! DO NOT allow children to play with packaging materials, as this presents a choking hazard if swallowed ! Carefully unpack your Deep Tissue Massager and check the contents for any missing items or any possible transport damage.

See "Product overview" section above for the contents that you should have.

If you do find any damaged items, then DO NOT use the product. Please return it to the place of purchase immediately with your receipt.

We recommend you retain the packaging for ease of storage and as a precaution in case you experience a problem and need to return the item for repair or exchange.

## CHARGE INSTRUCTIONS

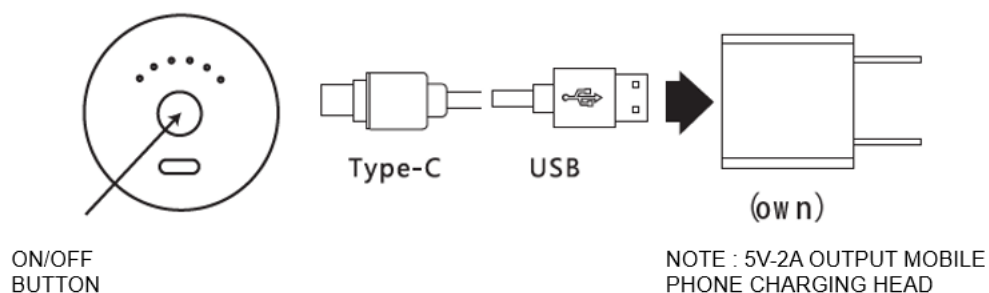1. Before the device is used for the first time or after not using it for a long period of time, it should be charged for 200 minutes (or till the 6 LEDs light-up permanently). 90 minutes cordless use time when battery is fully charged.
2. Long press for 3 seconds the ON/OFF button to switch off the massager prior to charging it with the provided USB charging cable.
3. Connect the USB charging cable into the device socket located on the bottom of the device's handle.
4. Insert the USB charging cable into a compatible adapter.
5. When charging the device: the 6 LEDs will alternatively flash to enable you to check the charging status.
6. When battery is fully charged: the 6 LEDs will light-up permanently.

7. Unplug the USB charging cable from the charging socket once the device is fully charged.

## OPERATING INSTRUCTIONS

### CORDLESS USE ONLY.

The massage treatment must feel pleasant at all times. A slight haematoma may well occur – this is a positive therapeutic effect as the use of the Deep Tissue Massager can selectively encourage your blood flow. In isolated cases, areas of redness may also occur on the skin.

**NOTE.** If you experience excessive skin irritation, stop using the massager and consult a doctor.

Depending on how tense/tight your muscles are, we recommend using the device for approximately 10-15 minutes on the area you wish to treat (legs, shoulders, arms and back).


## SELECTING / CHANGING THE MASSAGE HEAD

Before starting the massage, you first need to choose the right massage head. Section “PRODUCT OVERVIEW” contains a list of all the massage heads and what they are best suited for.

Long press for 3 seconds the ON/OFF button to switch off the device.

To change the massage head, proceed as follows:

1. PULL THE EXISTING MASSAGE HEAD HORIZONTALLY OFF THE DEVICE

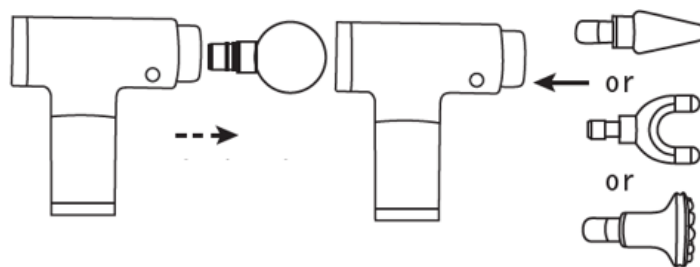2. INSERT THE DESIRED MASSAGE HEAD FIRMLY INTO THE DEVICE

Once the massage head is properly installed, long press for 3 seconds the ON/OFF button to switch on the device.

**GEAR ADJUSTMENT :** The starting gear is the 1<sup>st</sup> gear. To switch between the gears, briefly press the ON/OFF button once. This device features 6 gear settings.

**NOTE.** In order to protect the bones and muscles, avoid long-term fixed-point massage. massager comes with an automatic switch-off function. This means that the device automatically switches off after a period of use of 15 minutes.

## CARE AND CLEANING

**IMPORTANT:** ALWAYS ensure that the massager is switched off after use and allow it to cool down before cleaning.

**WARNING!** This device is NOT suitable for cleaning under an open water tap.

**WARNING!** DO NOT immerse the appliance into water or other liquids to clean it.

**IMPORTANT:** ALWAYS unplug the USB charging cable from the charging socket before cleaning.

DO NOT clean the appliance with steel wool, strong chemicals, alkali, abrasive or disinfecting agents as they may damage the surface.

Detach the massage head from the appliance if required.

DO NOT use the device again until it is completely dry.

## MAINTENANCE & TROUBLE SHOOTING

This appliance is designed for household use and requires minimal care after purchasing.

If the appliance malfunctions, switch it off and allow it to cool down, and contact your place of purchase for advice.

DO NOT attempt to repair the appliance yourself.

This appliance has no user-serviceable parts.

## STORAGE

When not in use, switch off the device. Allow the appliance to cool down completely and then store it out of the reach of children in a cool, dry location.

## ENVIRONMENTALLY FRIENDLY DISPOSAL

This symbol, displayed on the product or on its packaging, indicates that this product complies with the European Directive 2012/19/EU and must not be disposed of with household waste.

It must be taken to an appropriate collection point for recycling electrical and electronic equipment.

By ensuring that this product is disposed of in an appropriate way, you will be helping to prevent potential harmful consequences for the environment and human health.

Recycling materials will help to conserve natural resources.

Do not dispose of used batteries with household waste. They should be disposed of at a collection point. You can obtain more information from your local waste management company or your local council.

## WARRANTY

The warranty period starts from the date of purchase and is valid for a period of 24 months. The manufacturer is not liable for any damages other than those on products supplied.

### The following conditions are the warranty exceptions:

- The product is not used in accordance with the instruction manual.
- Users make alterations or modifications.
- Use after an accident.
- Damage due to unexpected disasters, such as fire, floods, etc.

**Official importer :**

**ZA, Les Petites Ruelles**

**28130 Saint-Piat**

**France**